

Specifications

Materials of Construction

Nasal Pillows: Medical grade silicone
 DreamSeal Mask: Perimeter seal is medical grade silicone
 Straps: Nylon and spandex
 All *Breeze SleepGear* components are latex free.

Pressure Range

3-35 cm H2O

Hose Attachment

Fits 22 mm OD connection



Breeze SleepGear With Nasal Pillows Assembly



Breeze SleepGear With DreamSeal Assembly



Nasal Pillows



DreamSeal Mask

Ordering Information

Description	Order Number
Breeze SleepGear Interface Systems With Nasal Pillows*	
Breeze SleepGear Interface System with Medium and Large Nasal Pillows	Y-101400-00
Breeze SleepGear Interface System with Small Nasal Pillows	Y-101400-S
Breeze SleepGear Interface System with Medium Nasal Pillows	Y-101400-M
Breeze SleepGear Interface System with Large Nasal Pillows	Y-101400-L
Breeze SleepGear Interface System with Extra Large Nasal Pillows	Y-101400-XL
Breeze SleepGear Interface System with Small Dilators	Y-101400-SD
Breeze SleepGear Interface System with Medium Dilators	Y-101400-MD
Breeze SleepGear Interface System with Large Dilators	Y-101400-LD
Breeze SleepGear Interface Systems With DreamSeal Mask*	
Breeze SleepGear Interface System with Standard DreamSeal Mask	Y-102616-00
Breeze SleepGear Interface System with Large DreamSeal Mask	Y-103059-00A
Breeze SleepGear Interface System with Shallow DreamSeal Mask	Y-103074-00A
Nasal Pillows Accessories/Replacement Parts	
Nasal Pillows Shell with pressure tap (for sleep lab use)	Y-101987-00
Nasal Pillows, Small, blue, 1 pair	S-616323-00A
Nasal Pillows, Medium, clear, 1 pair	S-616324-00C
Nasal Pillows, Large, green, 1 pair	S-616325-00A
Nasal Pillows, Extra Large, lavender, 1 pair	Y-101969-00
Dilator, Small, gray, 1 pair	S-616355-00C
Dilator, Medium, white, 1 pair	S-616356-00C
Dilator, Large, yellow, 1 pair	S-616357-00C
Nasal Airway Assembly without pillows, includes shell w/vent, 15-mm tubing, swivel adapter and hose guide	Y-102444-00
Nasal Pillows Shell, includes vent	Y-101699-00
Nasal Pillows Exhalation Vent, 3/pkg	Y-101974-00
Pillow Sizing Gauge, 10/pkg	S-616468-PG
DreamSeal Mask Accessories/Replacement Parts	
DreamSeal Mask Assembly, includes shell with vent, 15-mm tubing, swivel adapter, hose guide	Y-102620-00
DreamSeal Mask Assembly with seal, Standard	Y-102617-00
DreamSeal Mask Assembly with seal, Large	Y-103061-00A
DreamSeal Mask Assembly with seal, Shallow	Y-103076-00A
DreamSeal Shell, includes vent	Y-102622-00
DreamSeal Shell with pressure tap (for sleep lab use)	Y-102637-00
DreamSeal Exhalation Vent, 3/pkg	Y-102638-00
DreamSeal Mask, Standard	Y-102639-00
DreamSeal Mask, Large	Y-103060-00A
DreamSeal Mask, Shallow	Y-103075-00A
Breeze SleepGear Universal Components	
Breeze SleepGear Assembly, includes spine, cradle and outrigger	Y-102640-00
Breeze SleepGear Cradle Assembly	Y-102950-00
Halo Straps, 1 pair	Y-101975-00
15-mm Tubing and 15-mm/22-mm Swivel Adapter	Y-101976-00
Hose Guide Kit, includes hose guide, clip and halo straps	Y-102833-00
Swivel Adapter	S-616329-00B

* Interface Systems include halo straps, 15-mm tubing and airway assembly

This unique interface can be used with the Nasal Pillows or *DreamSeal™* Mask assembly.

Breeze SleepGear™

CPAP Interface System



Breathe easy.
 Sleep better.

Specifications are subject to change without notice.

Trademarks of Tyco Healthcare Group LP or an affiliate © 2004 Tyco Healthcare Group LP. All rights reserved.

C-AD-Breeze/GB

tyco

Healthcare

Puritan Bennett

Tyco Healthcare Deutschland GmbH
 Gewerbestraße 1
 D-93333 NEUSTADT/DONAU
 Tel.: +49 9445 9590
 Fax: +49 9445 959155

Tyco Healthcare Norden AB
 Gårdsvägen 18
 S-16967 SOLNA
 Tel: +46 (0) 8585 60 505
 Fax: +46 (0) 8585 60 506

Tyco Healthcare Spain S.L.
 Avda. San Pablo, 28
 E-28820 COSLADA, MADRID
 Tel.: +34 91 669 68 48
 Fax: +34 91 674 03 24

Tyco Healthcare Italia SpA
 Via Rivoltana, 2/d
 I-20090 SEGRATE (MI)
 Tel: +39 02 70317.1
 Fax: +39 02 7530125

Tyco Healthcare Nederland BV
 Hogeweg 105
 NL-5301 LL ZALTBOMMEL
 Tel: +31 (0) 418 57 66 00
 Fax: +31 (0) 418 57 67 93

Tyco Healthcare Belgium N.V.S.A.
 General de Wittelaan 9/5
 B-2800 MECHELEN
 Tel: +32 (0) 15 29 44 50
 Fax: +32 (0) 15 29 44 55

Tyco Healthcare Austria GmbH
 Campus 21, Europaring F09402
 A-2345 BRUNN/GEIRGE
 Tel: +43 (0) 2236 37 88 39
 Fax: +43 (0) 2236 37 88 40

Tyco Healthcare Switzerland Ltd.
 Roosstrasse 53, Postfach
 CH-8832 WOLLERAU
 Tel: +41 1 786 50 50
 Fax: +41 1 786 50 10

Tyco Healthcare France S.A.S.
 50, rue Pierre Curie
 78370 PLAISIR cedex
 Tel: +33 (0) 1 30 79 82 00
 Fax: +33 (0) 1 30 79 80 30

Tyco Healthcare Regional Export
 154 Fareham Road
 GOSPORT
 HAMPSHIRE
 UK - PO13 0AS
 Tel: +44 (0) 13 29 22 42 58
 Fax: +44 (0) 13 29 22 40 86

Tyco Healthcare Portugal
 Produtos de Saúde Ida.
 Est: Outeiro Polima, Lote 10-1º Piso
 Abóboda
 P-2785-521 S. Domingos de Rana
 Tel: +35 1 21 448 10 00
 Fax: +35 1 21 445 05 88

Tyco Healthcare South Africa
 Birchwood Court
 Montrose Road, P.O. Box 8108
 Halfway House
 1685 MIDRAND
 SOUTH AFRICA
 Tel: +27 (0)11 265 2000
 Fax: +27 (0)11 265 2050

Tyco Healthcare South Africa
 Birchwood Court
 Montrose Road, P.O. Box 8108
 Halfway House
 1685 MIDRAND
 SOUTH AFRICA
 Tel: +27 (0)11 265 2000
 Fax: +27 (0)11 265 2050

Tyco Healthcare UK Limited
 154 Fareham Road
 GOSPORT
 HAMPSHIRE
 UK - PO13 0AS
 Tel: +44 (0) 13 29 22 41 59
 Fax: +44 (0) 13 29 22 40 83

07/2004

tyco

Healthcare

Puritan Bennett

Breeze SleepGear™ Interface

Your First Choice in Fit and Comfort.


Breeze SleepGear CPAP Interface System

The unique *Breeze SleepGear*™ CPAP Interface System from Puritan Bennett sets new standards in comfort and versatility for sleep apnea patients. Unlike traditional headgear and masks, there are no complicated or confining straps to contend with. The streamlined design contours the head and floats over the face, giving patients the freedom to sleep in virtually any position.

Two-in-One Versatility

The innovative *Breeze SleepGear* system lets patients easily switch between the Puritan Bennett Nasal Pillows and *DreamSeal*™ Mask, so they aren't limited to the same type of interface night after night. The patented Nasal Pillows are available in seven sizes, offering a wide choice in fit. The *DreamSeal* Mask comes in three sizes.

For more information visit www.puritanbennett.com or www.tycohealthcare.com

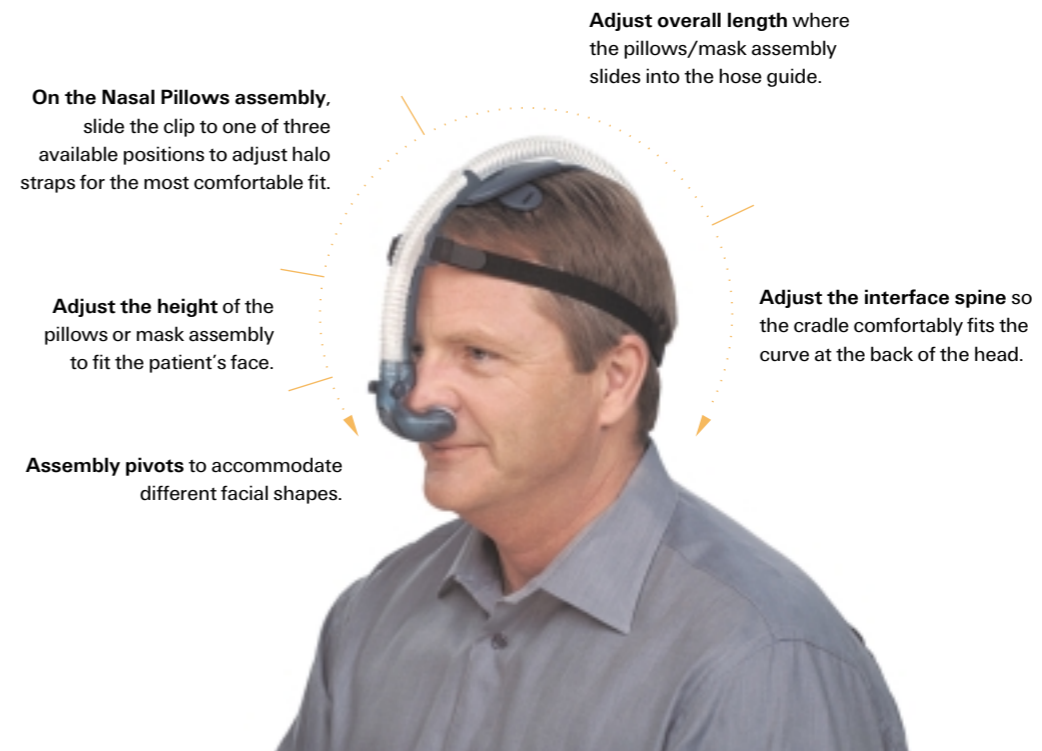

- > New clip (Nasal Pillows only) for patients who want more stability.
- > Multiple adjustments and size options make it easy to get just the right fit.
- > Versatile headgear works with the Nasal Pillows or *DreamSeal* Mask.
- > No confining, hard-to-adjust straps.
- > Ultra-quiet exhalation vent.
- > Patients can sleep in virtually any position.

New & Improved Features:

For over five years, users of the *Breeze SleepGear* CPAP Interface System have shared their gratitude for its comfort and performance. They also shared feedback and suggestions that we listened to in adding new features to make *Breeze SleepGear* system even better.

- > **Durability**—The hose guide for the Nasal Pillows assembly has an improved design and features a stronger material to ensure greater longevity.
- > **Fit**—A new clip has been added to the Nasal Pillows assembly to give patients greater flexibility in positioning the halo straps for a secure fit.
- > **Comfort**—The cradle pad has thicker padding to increase comfort for a good night's sleep.
- > **Convenience**—A cuff has been added to the swivel adapter for improved convenience.

Fitting the Breeze SleepGear Interface with Nasal Pillows



Make the *Breeze SleepGear* system your first choice in headgear for its ability to accommodate a wide range of patients. Simply slide the Nasal Pillows or *DreamSeal* Mask assembly into the *Breeze SleepGear* interface, and take advantage of the key adjustment points shown below.